




REVERE[®] ADDITION[®]

Revision Fixation System

Revision Surgery

Multiple Challenges, One System

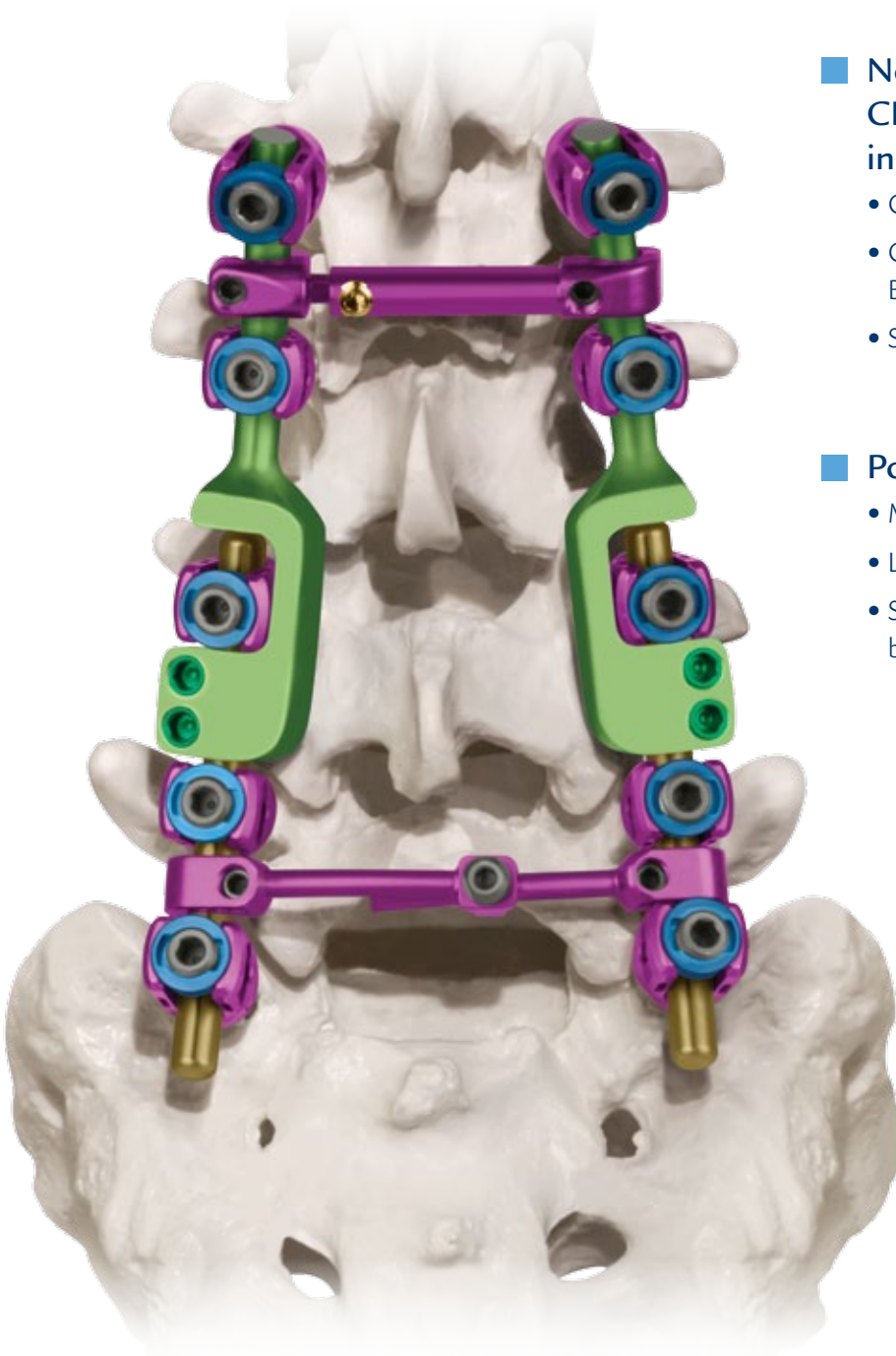

■ No need to remove existing REVERE®, CREO®, BEACON® or PROTEX® instrumentation


- Comprehensive range of implants
- Compatible with Globus' REVERE®, CREO®, BEACON® and PROTEX® systems
- Solid connection at critical junction

■ Potentially lower patient morbidity

- Multiple revision options
- Less scar tissue exposure
- Smaller incision may help to reduce blood loss and overall surgery time

Wide Variety of Unique Connectors

RESCUE® Implants

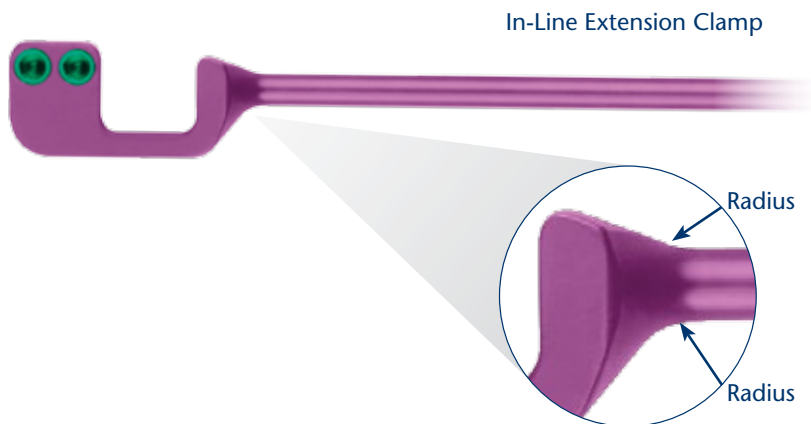

In-Line Extension Clamp

In-Line Extension Clamp

- Attach to existing Globus' REVERE®, CREO® BEACON®, or PROTEX® hardware
- Provide direct in-line coupling
- Integrated rod for stabilization
- Easier access to tighter areas (ideal for smaller stature patients)

The Strength is in the Radius

The critical junction where the rod meets the clamp (the radius) is mechanically stronger than a typical revision connector.



Offset Extension Clamp

Offset Extension Clamps

- Attaches to existing REVERE®, CREO®, BEACON® or PROTEX® hardware
- Provides 7mm offset coupling
- Integrated rod for extension



"U" Extension Clamp

"U" Extension Clamps


- Shorter run-on-rod with built-in lordotic curve
- Especially good for lumbar connections
- Offsets medially or laterally depending upon anatomical requirements
- Integrated rod

"Y" Clamps

- Shorter run-on-rod on each side of the screw
- Offsets medially or laterally, to meet anatomical requirements



"Y" Clamps

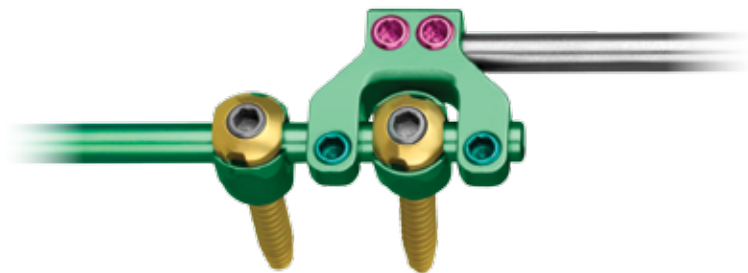

Other Available Implants



Parallel Connector Clamps




Rod-to-Rod Clamps

Utilizing the Unique REVERE® Non-Threaded Locking Cap


- Non-threaded design eliminates cross-threading
- 90° rotation of locking cap captures the rod
- Set screw assembly allows for easy insertion and tightening



Locking Cap

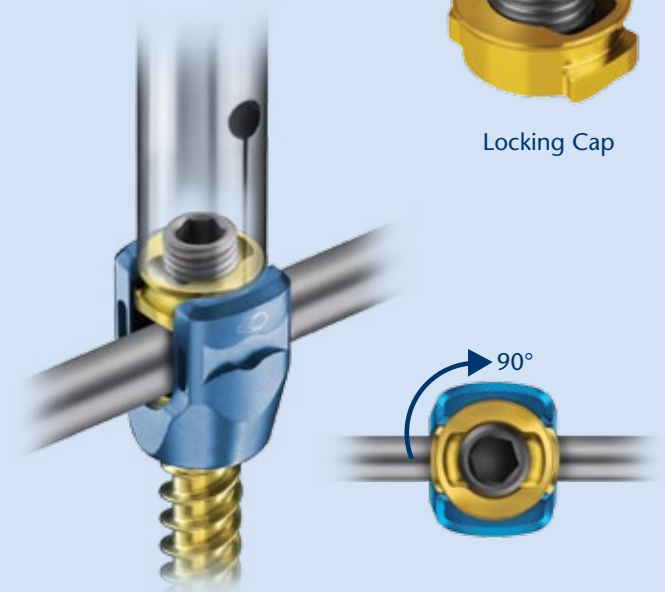

Locking Cap
INSERTED

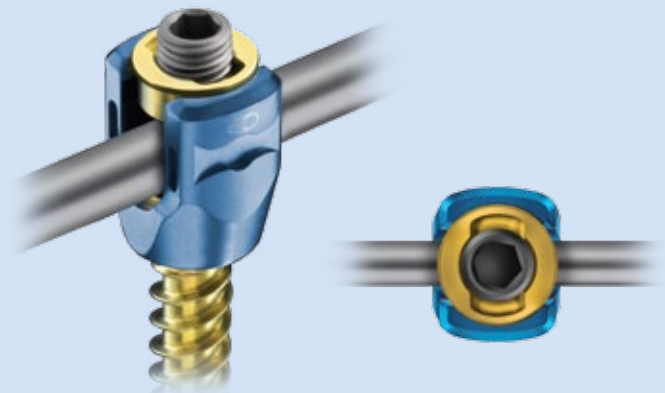

Locking Cap
LOCKED


Surgeon Testimonials

- *"I have never seen anything like this, I have another revision this Thursday, can you bring it...I will switch the case."*
—Dr. G, Detroit, MI
- *"With the most comprehensive selection of configurations and sizes, the Globus REVERE® ADDITION® Revision System has proven to be invaluable in helping me with my most challenging revision patients."*
—Dr. C, Loma Linda, CA
- *"I appreciate the intraoperative options that the REVERE® ADDITION® set offered me and that I didn't have to perform a longer exposure."*
—Dr. C, Dallas, TX



REVERE® ADDITION® Advantages

- Fits 5.5–6.5mm rods
- Lower overall implant costs
- Ease of use

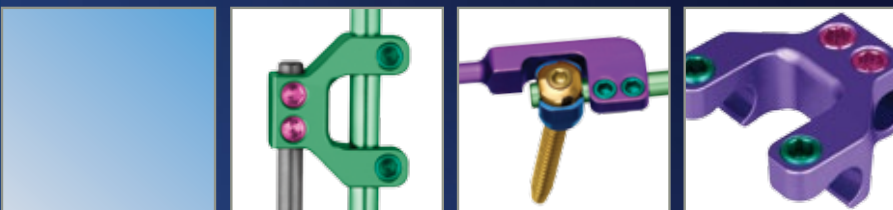

Globus Medical
Valley Forge Business Center
2560 General Armistead Avenue
Audubon, PA 19403
www.globusmedical.com

Customer Service:

Phone 1-866-GLOBUS1 (or 1-866-456-2871)

Fax 1-866-GLOBUS3 (or 1-866-456-2873)